



Radiation Monitors

www.lablogic.com

 **LabLogic**
EXPERIENCE & EXPERTISE

Radhound Monitors

Radhound monitors are ideal for use within a large range of sectors, including academia, pharmaceutical research and nuclear medicine/PET. These meters are compatible with a broad range of Geiger-Müller, sodium iodide and plastic scintillator based probes to cover the whole spectrum of contamination and dose rate measurement requirements.

Radhound Monitors	3 - 4
Contamination Probes	5
Dose Rate Probes	6
Scintillation Probes	7
Alpha / Beta Contamination Probes	8
Handhound	9

Tracerco Monitors

The Tracerco range is unique, robust and innovative – great for use in the laboratory or by radiation protection professionals out in the field.

Contamination, Dose Rate and X-ray Monitors	10
Personal Electric Dosimeter (PED)	11
NORM Monitor	11

Radhound X/E and X/I

The Radhound X/E is an advanced hand-held general purpose radiation monitor, suitable for a wide range of probes. The X/I is a Radhound with an internal dose rate detector.

This feature-packed instrument boasts some unique features, such as the ability to switch between probes via the menu allowing, for example, a dose rate probe and a contamination probe to be configured for use with one instrument. This flexibility allows any standard probe to be used (operating voltage 300 - 1200 V).

- Clear digital LCD display with backlight.
- GM and scintillation detector options.
- Multiple probe library/configuration.
- Fully adjustable alarm levels.
- Scaler timer function.
- Over range indication.
- Peak mode.

Radhound Multi-purpose Digital Radiation Meter

A multi-purpose digital radiation survey meter suitable for all your contamination monitoring and radiation protection requirements, the Radhound is a cost effective, feature packed digital radiation monitor that is simple and easy to use.

Count rate is displayed in large clear numbers and also on a bar scale. Our smart averaging software means a steady display that can be read with confidence, yet provides a fast response.

For source finding, one button push changes the display to a histogram plot. Alpha and Beta/Gamma counts can be displayed separately or on the same screen.

For surveying operations the Radhound also has an integrator mode.

- Clear digital LCD display with backlight.
- GM and scintillation detector options.
- Scaler timer function.
- Ergonomic tilt stand.
- Wall mountable.
- Fully adjustable alarm levels.



Radhound X335 Advanced Hand-held Survey Meter

The Radhound X335 is an advanced hand-held survey meter with an integrated compensated end window Geiger-Müller detector for gamma dose rate monitoring.

- Dose Rate Range H*(10): 0 uSv/h – 1 mSv/h
- Energy response 20 keV - 1.5 MeV
- Clear digital LCD display with backlight.
- Adjustable alarm levels.
- Timed count and timed dose function.
- Peak mode.
- Over range indication.



Radhound Alarm

A low cost digital area monitor available with a range of probe options for all your contamination monitoring and radiation protection requirements.

The Radhound Alarm has a large digital display, the option for mains or battery operation and can be easily wall mounted. The easy to use password protected configuration menu enables a choice of probe and alarm thresholds to be setup.

Alarm output is provided by remotely mounted audio sounder and visual beacon.

- Simple to operate.
- Programmable alarm threshold.
- Audible count rate (with on/off switch).
- Customisation available.



SS300 Probe

The SS300 is an uncompensated pancake Geiger-Müller-based probe for alpha beta and gamma contamination measurement.



SS315 Probe

The SS315 is functionally identical to the SS300, but with a different probe geometry.



Efficiencies

(Listed as percentage of 2π emission rate)

Nuclide	Emission	Efficiency
²⁴¹ Am	α	29.1%
²³⁸ Pu	α	26.6%
Nat U	α	63.5%
⁹⁰ Sr/ ⁹⁰ Y	β	56.7%
¹⁴ C	β	19.4%
¹⁴⁷ Pm	β	59.1%
²³⁸ Pu	β	25.8%
⁶⁰ Co	β	36.2%
¹³⁷ Cs	β	50.6%

	SS300 (Pancake)	SS315 (End On)
Operating Voltage	550 V	550 V
Measurement Range	0 - 5 kcps	0 - 5 kcps
Plateau Length	150 V	150 V
Temperature Range	-10°C to + 50°C	-10°C to + 50°C
Energy Response	20 keV - 2 MeV Gamma, ≥ 40 keV Beta, > 3 MeV Alpha	20 keV - 2 MeV Gamma, ≥ 40 keV Beta, > 3 MeV Alpha
Connector Type	MHV	MHV
Dimensions	∅ 70 x 254 x 64 mm	∅ 70 x 180 mm
Active Area	15.5 cm ²	15.5 cm ²
Mass	280 g	450 g

Dose Rate Probes

SS330 Probe

The SS330 probe is an excellent general purpose end window compensated pancake Geiger-Müller probe with H*(10) energy compensation, which permits reliable measurements from ambient background up to 1 µSv/hr.

SS335 Probe

The SS335 probe is functionally identical to the SS330, but with a different probe geometry.

SS340 Probe

The SS340 is a side-window Geiger-Müller probe for ambient gamma radiation measurement to H*(10).

Dose rate range is 0 - 2 mSv/hr and energy range 45 keV - 2 MeV.



	SS330	SS335	SS340
Operating Voltage	550 V	550 V	450 V
Dose Rate Range*	0.1 µSv/hr - 1 mSv/hr	0.1 µSv/hr - 1 mSv/hr	0.1 µSv/hr - 1 mSv/hr
Plateau Length	150 V minimum	150 V minimum	200 V minimum
Temperature Range	-10°C to + 50°C	-10°C to + 50°C	-10°C to + 50°C
Gamma Sensitivity	Typically 5 cps/µSv/hr	Typically 5 cps/µSv/hr	Typically 2 cps/µSv/hr
Energy Range	H*(10) for 20 keV - 1.5 MeV	H*(10) for 20 keV - 1.5 MeV	H*(10) for 45 keV - 1.5 MeV
Connector Type	MHV	MHV	MHV
Dimensions	Ø 70 x 254 x 64 mm	Ø 70 x 180 mm	Ø 25 x 135 mm
Active Area	15.5 cm ²	15.5 cm ²	40 mm tube length
Mass	300 g	470 g	100 g

* Dose rate probes are set up to read in µSv/hr by default. For measurements in rem/hr, please specify at point of order.

Scintillation Probes

SS404 Al Probe and Be Probe

The SS404 Al is a thin-crystal NaI(Tl) end-window scintillation probe designed to be an equivalent to the Mini 44A.

This probe incorporates a Ø 1.25" x 0.098" (32 x 2.5 mm) thick NaI(Tl) crystal mounted on an aluminium window and is fitted with an internal 3.15 mm lead collimator to reduce background counts.

The SS404 Be is similar to the SS404 Al but is fitted with a beryllium window, which extends the low energy response down to 5 keV, making it suitable for counting ⁵⁵Fe.



SS500 Probe



The SS500 is a very sensitive end-window gamma scintillation probe.

Equipped with a Ø 1" x 1" (25.4 x 25.4 mm) NaI(Tl) crystal, it is designed to provide a cost effective gamma monitor for energies of 50 keV upwards.

	SS404 Al	SS404 Be	SS500
Operating Voltage	Typically 650 V	Typically 650 V	Typically 650 V
Detector Crystal	Ø 1.25" x 0.098" NaI (32 x 2.5 mm)	Ø 1.25" x 1" NaI (32 x 2.5 mm)	Ø 1" x 1" NaI (25.4 x 25.4 mm)
Window Density	2.7 mg/cm ³	1.85 mg/cm ³	35 mg/cm ²
Energy Response	15 keV - 250 keV	5 keV - 250 keV	50 keV - 2.0 MeV
Connector Type	MHV	MHV	MHV
Dimensions	Ø 54 x 185 mm	Ø 54 x 185 mm	Ø 44.5 x 205 mm
Mass	820 g	820 g	300 g
Operating Temperature	-10°C to + 50°C	-10°C to + 50°C	-10°C to + 50°C
Operating Humidity	Up to 95% RH non-condensing	Up to 95% RH non-condensing	Up to 95% RH non-condensing

Efficiencies (Listed as percentage of 2π emission rate)

Nuclide	Energy	SS404 Al Efficiency	SS404 Be Efficiency
⁵⁵ Fe	5.9 keV	-	31.4%
²³⁸ Pu	16.3 keV	98.7%	99.1%
¹²⁹ I	31.5 keV	84.9%	91.5%
²⁴¹ Am	58.8 keV	117.0%	117.3%
⁵⁷ Co	120 keV	82.7%	83.0%
¹³⁷ Cs	662 keV	17.0%	18.3%
⁶⁰ Co	1200 keV	11.4%	12.4%

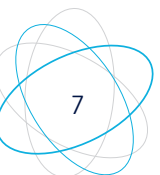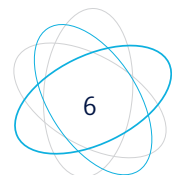


SS600 Probes

Equivalent to the NE BP6 / AP2, there are three versions of these 100 cm² window probes available:

- SS600 A Alpha only (Zinc sulphide layer)
- SS600 B Beta only (Plastic scintillator)
- SS600 AB Alpha/Beta (Zinc sulphide bonded to a plastic scintillator)

The use of a plastic scintillator avoids the traditional use of anthracene in this application, with a comparable response.

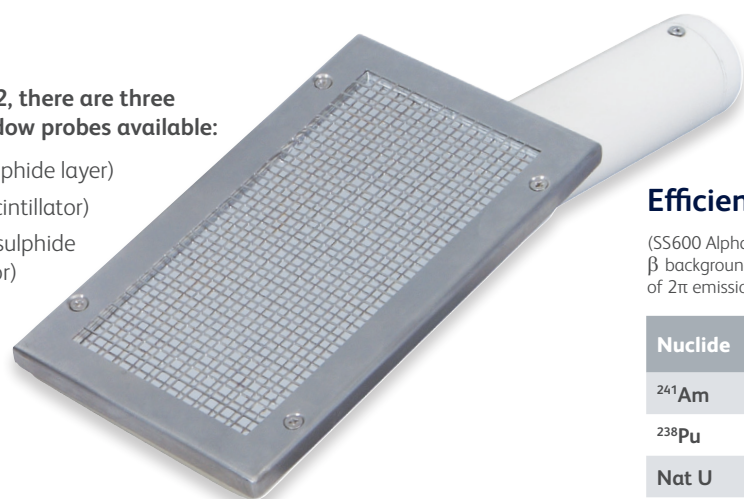

SS700 Probes

A series of three ergonomically balanced probes with a square window of 50 cm² and a 64° angled handle.

Equivalent to the NE BP7, there are three versions available:

- SS700 A Alpha only (Zinc sulphide layer)
- SS700 B Beta only (Plastic scintillator)
- SS700 AB Alpha/Beta (Zinc sulphide bonded to a plastic scintillator)

The use of a plastic scintillator avoids the traditional use of anthracene in this application, with a comparable response.



Efficiencies

(SS600 Alpha/Beta, α background 1.9 cps, β background 7.5 cps. Listed as percentage of 2 π emission rate)

Nuclide	Emission	Efficiency
²⁴¹ Am	α	39.3%
²³⁸ Pu	α	42.0%
Nat U	α	43.1%
⁹⁰ Sr/ ⁹⁰ Y	β	38.4%
¹⁴ C	β	1.5%
³⁶ Cl	β	36.9%
²³⁸ Pu	β	4.7%
⁶⁰ Co	β	14.0%
¹³⁷ Cs	β	28.8%

Efficiencies

(SS700 Alpha/Beta, α background 0.7 cps, β background 3.4 cps. Listed as percentage of 2 π emission rate)

Nuclide	Emission	Efficiency
²⁴¹ Am	α	33.0%
²³⁸ Pu	α	32.1%
Nat U	α	34.8%
⁹⁰ Sr/ ⁹⁰ Y	β	33.9%
³⁶ Cl	β	30.6%
⁶⁰ Co	β	13.2%
¹³⁷ Cs	β	23.1%

Alpha

	SS600 A	SS700 A
Window Area	100 cm ²	50 cm ²
Typical Alpha Efficiency ²⁴¹ Am % of surface emission (2 π)	30	30
Background Response	Less than 1 cps	Less than 1 cps
Mass	500 g	500 g

Beta

	SS600 B	SS700 B
Window Area	100 cm ²	50 cm ²
Beta Efficiency ¹⁴ C % of surface emission (2 π)	8	40
Beta Efficiency ⁹⁰ Sr/Y % of surface emission (2 π)	8	40
Background Response	Less than 20 cps	Less than 10 cps
Mass	500 g	500 g

Alpha / Beta

	SS600 AB	SS700 AB
Window Area	100 cm ²	50 cm ²
Typical Efficiency ²⁴¹ Am % of surface emission (2 π)	30	30
Typical Efficiency ⁶⁰ Co % of surface emission (2 π)	15	15
Typical Efficiency ⁹⁰ Sr/Y % of surface emission (2 π)	30	30
Background Response	Less than 20 cps	Less than 10 cps
¹³⁷ Cs Response cps/ μ Sv/hr	50	25
Mass	500 g	500 g

Handhound Voice Activated Monitor

The Handhound has been designed mainly for use in 'wet chemistry' radio-isotope handling situations where hands could be contaminated. It has been designed to help with HSE compliance, by keeping a record of hand contamination measurements taken.


Radhound radiation monitoring electronics are incorporated into a stainless steel enclosure with a high sensitivity collimated detector, so that any contamination can be detected on hands.

The background is updated while the unit is not in use. A proximity sensor ensures the user's hands are underneath the detector, the user then speaks his/her name and says 'Continue'. The system will begin counting for a predetermined period, set by the supervisor.

A touchscreen interface is also incorporated, to allow configuration and manual triggering if needed.

The Handhound is mains-operated and will find application in 'wet' lab using gamma isotopes, in corridors, work areas and laboratories.

- Entirely voice operated to avoid instrument contamination.
- Sensitive scintillation counter for gamma emitters.
- Alternative detector options covering a wide range of nuclides.
- Stainless steel housing for ease of cleaning and decontamination.
- Automatic record keeping against user names, to aid with HSE compliance.
- Touch-screen compatibility included as an alternative to voice operation.
- Data can be downloaded onto USB.
- Automatic background updates.
- Fixed or dynamic alarm thresholds.
- Mains operated.

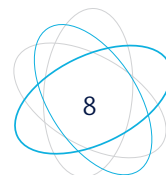



T402 and T406 Dose Rate and X-ray Monitors

The T402 and T406 are lightweight, yet robust and comfortable to use over extended periods.

- T402 detects gamma and X-rays from 60 keV - 1.33 MeV.
- T406 detects gamma and X-rays from 17 keV - 1.33 MeV.
- T402HR – extended range for high dose rates.
- Digital bar graph display: 0.1 - 1000 $\mu\text{Sv/h}$.
- Digital dose rate indication: 0 - 10,000 $\mu\text{Sv/h}$.
- Peak dose rate memory – allows maximum exposure levels to be recorded.
- Accumulated dose rate memory – for risk assessment and total exposure.
- Audible response with adjustable alarm thresholds.
- Water-resistant so easy to clean and decontaminate.
- Shock and drop tested so highly durable.



T401 Contamination Monitor

Designed to meet the challenge of combining operational reliability with excellent sensitivity the T401 offers a range of features including direct surface, peak and background readings. It can be used one-handed, or detach the probe for two-handed operation.


The T401 can be supplied with an extension pole kit to securely deploy the detector probe during monitoring operations.

- Dual bar graph meter display 0 - 1000 cps.
- Digital numeric display with automatic direct translation to Bq/cm^2 for 14+ pre-programmed nuclides (natural and man-made) including ^{14}C , ^{32}P , ^{137}Cs .
- Optional extension arm.
- Detachable probe.
- Background reading and storage.
- Audible response with adjustable alarm thresholds.

Personal Electronic Dosimeter (PED)

Ideal for users who are not specially trained to measure radiation exposure, the PED family have been specially designed to be easy to use and understand. Encased in weather, shock and drop-proof housings each PED features a smooth clean design and simple to use DoseVision™ software.

- Detects X-rays and gamma rays from 48 KeV - 3MeV.
- One button operation.
- Easy to read large AMOLED display screen displaying dose rate, accumulated dose and animated silhouette indicating dose received.
- Multiple languages.
- Multiple users.
- Waterproof up to 1 m.

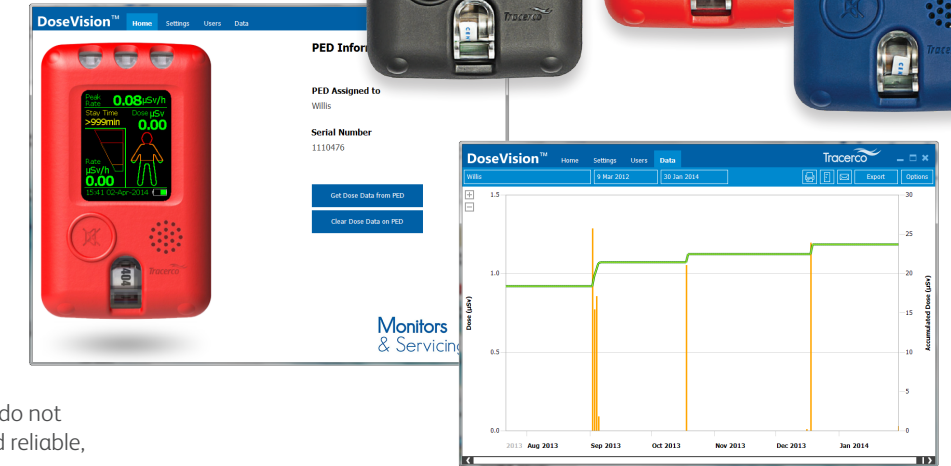

PED-IS

This intrinsically safe PED is perfect for both radiation specialists and those who do not work with radiation every day. Robust and reliable, it is safe to use in potentially explosive areas, making it ideal for challenging environments.

PED-Blue

This is the non-intrinsically safe version of the PED-IS. Lighter, it retains the same high quality design and features a direct micro USB connection.

This entry level model is perfect for use within the medical sector.

PED+

An advanced version of the PED-Blue, it can be used as both a PED and a hand-held dose rate survey meter. The PED+ has a number of added features, such as Bluetooth, GPS and pop-up message alarms.

PED-ER

The PED-ER's extended dose rate range of 1 Sv/h (100 R/h) provides perfect radiation dosimetry for nuclear medicine environments.



PED ER+

Featuring an extended dose rate range of 1 Sv/h (100 R/h), the PED-ER+ allows the user to measure radioisotopes which could not have been measured previously. Lightweight, waterproof and compact, the PED-ER+ provides the perfect radiation monitoring solution for those working in challenging environments.

NORM Monitor-IS

The ultimate tool for obtaining accurate NORM measurements in hazardous areas, the NORM Monitor-IS is ATEX approved with dual probe capabilities; Geiger-Müller and scintillator.

- Large, easy to read LCD screen with bar graph and back light.
- One-touch integrate function that allows detection of very low activities for increased measurement accuracy.
- Live background subtraction and several measurement modes.
- Adjustable alarm thresholds for improved safety.
- Easy to clean and decontaminate.

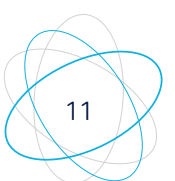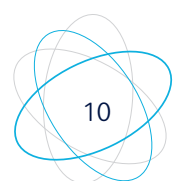

Service and Support

Users of our systems can benefit from our comprehensive, fully inclusive service and support.

We can give reassurance that if things go wrong or you need expert advice, help is only an e-mail or phone call away.

Validation Services

Our Validation Service enables you to implement and get maximum value from your investments as soon as possible.

We work as a partner with your Quality Manager, System Manager and users to provide a tailored Validation Plan, suited to your needs. Our Validation Specialists have years of experience in GLP system validation, detailed knowledge of our systems, together with other industry standard systems to help you meet company and regulatory requirements.

Training

LabLogic can provide a variety of training courses and workshops to help you get the most out of your instrument and software.

All training is performed by our expert Product and Support Specialists who have many years experience in the development and use of the instruments and software.

Certificates can be provided to complement your internal GLP training records.

Europe & Worldwide

LabLogic Systems Limited

Paradigm House, 3 Melbourne Avenue
Broomhill, Sheffield, S10 2QJ, UK

E-mail: solutions@lablogic.com

Tel: +44 (0)114 266 7267

Fax: +44 (0)114 266 3944

www.lablogic.com

USA & Canada

LabLogic Systems, Inc.

1911 N US HWY 301, Suite 140
Tampa, FL 33619, USA

E-mail: solutions@lablogic.com

Tel: +1-813-626-6848

Fax: +1-813-620-3708

www.lablogic.com



LabLogic

EXPERIENCE & EXPERTISE



Certificate No: 1535
ISO 9001



INVESTORS
IN PEOPLE

Silver
Until 2022



Certificate No: 10926
ISO 9001