



RADIOPROTECTION:
IT IS THERE FOR YOU AND IT IS EVIDENT.

KARL100 is the smallest and safest Class A radiopharmaceutical dose administration system available on the international market, with 4 major benefits.





# **OPERATOR SAFETY**

NO HAZARDOUS MANUAL OPERATIONS



Uncapping of the lead pot and perforation of the multi dose vial septum are performed automatically, as well as the dose drawing from the vial, within the fully shielded compartment for a full operator protection throughout the whole process. A specific study performed with an innovative dose rate monitor (Aegis ED3) proved a significant decrease of the staff doses using KARL100.

### DOSE MAP

The highest dose received to each part of the hand is rendered in red, and the lowest in blue.





**KARL100** is also equipped with a built-in GM tube for a **real-time online dosimetry monitoring**. **Remote control** of the system via any commercial tablet is also (optionally) available.





## PATIENT SAFETY

### KEEPS RADIOPHARMACEUTICALS STERILE



Your patient is never connected to the multi-dose mother vial: there is **no risk of patient-to-patient contamination**. Each single syringe is automatically connected to the main disposable set within an isolated vented box with a filtered and controlled atmosphere (Class A). Incoming activity from multi-dose vial is filtered through an **inline 0.22 micron filter**.

#### NO RISK OF AN OVERDOSE

The injection of single dose syringes by means of **mini injectors** prevents any chance of an accidental dangerous overdose.






# **FLEXIBILITY**

### USE YOUR STANDARD SHIELDED CONTAINERS



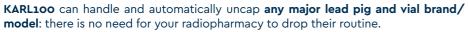

#### NO LIMITS FOR SUPPLIED ACTIVITY CONCENTRATIONS



**KARL100 draws all of the source vial activity up to the last drop!** Also, it grants a clear-cut dispensing, even for an extremely high incoming activity concentration, thanks to an automatic dilution of the main activity.





#### CARRY ONLY 1KG MINI-INJECTORS!



Quit carrying heavy shields each time you inject a dose: the operator may carry only light and portable mini-injectors to each injection room (user-friendly), allowing dose preparations and multiple infusions simultaneously (time-friendly) and reducing patient anxiety during dose administration by using only these very compact individual mini-injectors (patient-friendly). It is always possible to measure a dispensed dose into an external calibrator (ie. for quality control) without having to waste it. The dose remains sterile and can be injected right afterwards.











#### THE PET INFUSION REVOLUTION

DISCOVER ALL THE BENEFITS OF THE SMALLEST AND SAFEST AUTOMATED ADMINISTRATION SYSTEM.

#### **DIMENSIONS & CHARACTERISTICS**

#### Mini Injector

Weight: 1.1 Kg (2.4 lb) Length: 16.3 cm (6.4")

Class 2b (CE 0476) Medical Device Power Supply: Lithium Battery

#### **Dose Drawing Machine**

Weight: 445 Kg (981 lbs)

Dimensions: 54(w) x 116(d) x 107(h) cm (21" x 45" x 42") Power Supply: 110/220 VAC 50/60 Hz + Built-in Lithium Battery

Power Requirement: Max 300 Watts Class 1m (CE 0476) Medical Device

#### KARL100 FEATURES

- Automatic connection of the syringe to the main disposable set within a GMP Class A, HEPA filtered atmosphere
- Incoming activity from multi dose vial filtered through an in-line 0.22 micron filter
- No risk of patient-to-patient contamination / No risk of accidental overdose
- Handling and automatic uncapping of the lead pig for any major brand pig/vial model combination
- Automatic perforation of the multi dose vial septum

DOSIMETRY DATA AT REST (in μSv/h):				
Incoming PET activity	Surface	5 cms from Surface	10 cms from Surface	30 cms from Surface
700 mCi	12,00	7,11	4,72	1,0
750 mCi	12,78	7,62	5,06	1,07
800 mCi	13,63	8,13	5,40	1,14
1 Ci	16,00	8,5	6,47	1,5
1,5 Ci	25,5	15,20	10,12	3,33
2,0 Ci	34,70	20,03	13,49	4,44

- Max concentration of incoming activity: no limits
- Max incoming PET Radiopharmaceutical volume: 40 mL
- Volumetric Dose Dispensing Precision: +/- 2.0 %
- Provided with built-in GM tube, for continuous record of dosimetry safety monitoring GAMP5, 21 CFR part11 compliant control software provided with:
  - User-Friendly touch-screen total control interface
- Configurable weight based dose calculation algorithm for patient dose optimization
- Cabled / WiFi Ethernet connection for allowing efficient internet communication for trouble-shooting and (Optional) data exchange with PET/NM departmental management software
- Remote control of the system, by use of a (Optional) commercial tablet
- Motorized Drive for moving the dose drawing machine easily
- Battery with 8-hours working time with 100% recharging in only 2 hours!







# MULTI-TRACER COMPATIBLE

<sup>18</sup>F - <sup>11</sup>C - <sup>177</sup>Lu - <sup>90</sup>Y <sup>223</sup>Ra - <sup>68</sup>Ga - <sup>64</sup>Cu



