

Hidex AMG

Automatic Gamma Counter



HIDEX

LabLogic

EXPERIENCE & EXPERTISE

Effortless operation with integrated applications

The Hidex AMG is the first gamma counter dedicated to Nuclear Medicine and PET applications.

Specifically designed to meet the needs of modern laboratories, the Hidex AMG's touchscreen operated software and application-focused design delivers effortless workflow and provides results quickly and efficiently.

Facilitating efficiency and workflow

- Combine all the necessary individual measurement data in a glomerular filtration rate (GFR) assay to produce a final report for secure workflow in hospitals.
- Add easy quality control through Ge/Ga elution ratio measurements.
- High efficiency for PET research with isotopes such as ^{15}O , ^{18}F and ^{89}Zr with automatic decay correction to a defined time stamp.



- An easy to use diagnostic tool for traditional radio immunoassays with automatic data export into user definable format.
- Ensure radio safety through ready-set gamma counting measurements.


Technical Innovations

A unique combination with the onboard balance

The optional automated onboard balance by Sartorius, offers unprecedented convenience and reliability.



It saves the operator valuable time by speeding up work flow. Results are reported automatically as activity per mass or volume providing increased dependability.

The Hidex AMG with the precise and easy to use onboard balance is a perfect platform for bio distribution studies.

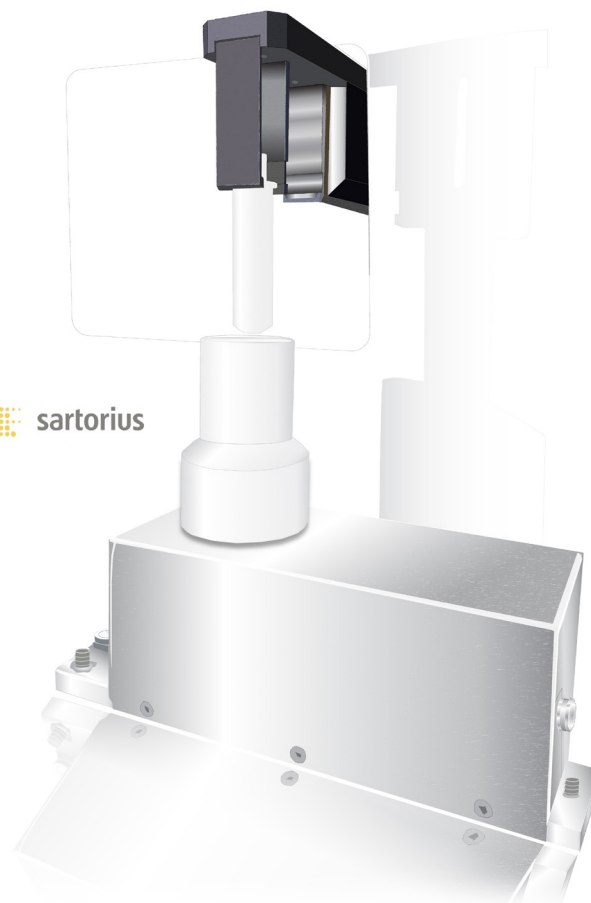

Going beyond with the extended MCA

The Hidex AMG is equipped with a powerful linear multichannel analyzer for detailed spectrum analysis. Never before has there been a possibility to measure isotopes with an energy range of up to 4000 keV.

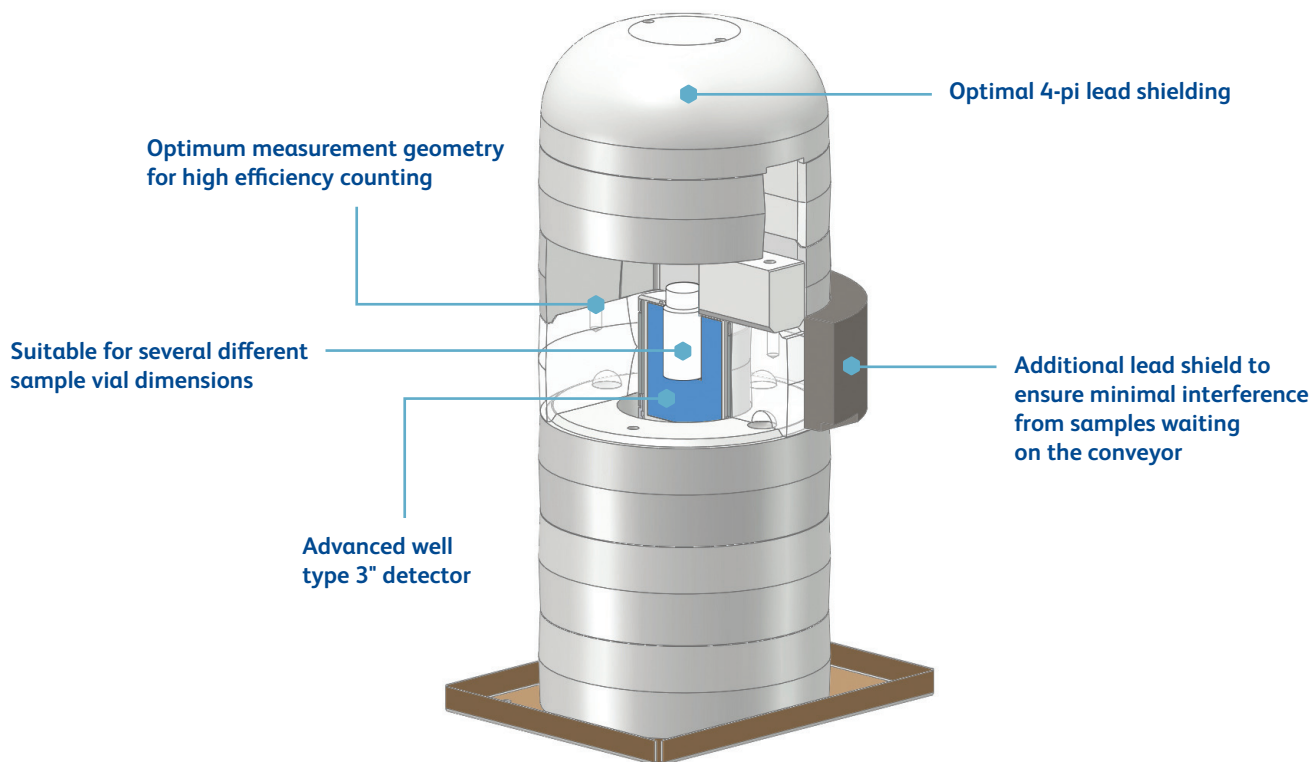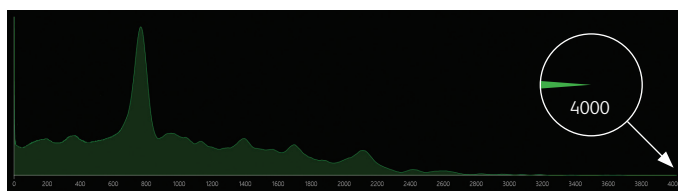

Software beyond comparison

The Hidex AMG is ready to run with immediate access to your applications. Results are calculated straight into Bq or Ci values and automatic storage of the full energy spectrum ensures no data is lost.

An external PC with Windows 10 operated software features an easy to use graphical interface, with data exporting to a user defined location.

The software has an in-built algorithm to compensate for spectral pulse pile up of high activity samples.

Besides traditional printing, it is also possible to store electronic data in several formats. Sample ID input ensures easy handling of data.

In order to satisfy laboratory security demands, the Hidex AMG software is 21 CFR part 11 supported with user rights and permissions combined to Windows user accounts.



Specifications

Size	25" (W) x 27" (D) x 23" (H) (620 x 700 x 600 mm)
Weight	451 lb (198 kg)
Sample Quantity	Ø 13 mm max 250 pcs. and Ø 28 mm max 78 pcs.
Counting Efficiency (Peak Window)	Cs-137 18%, I-129 58%
Background (Open Window, CPM)	350

Please refer to the Technical Specification Sheet for further information



Specifications

Size	25" (W) x 28" (D) x 24" (H) (620 x 700 x 600 mm)
Size with balance	25" (W) x 33" (D) x 24" (H) (620 x 830 x 600 mm)
Weight	437 lb (198 kg)
Weight with balance	448 lb (203 kg)
Sample Capacity	Ø 13 mm max. 250 pcs. and Ø 28 mm max. 78 pcs.; 95 mm max. height
Counting Efficiency	Cs-137 47%, I-129 58%, Cr-51 6%
Detector Shielding	Optimal 4-pi Pb shield design 55 mm; 80 mm shielding on conveyor side
Background (Open Window)	< 328 cpm

Please refer to the Technical Specification Sheet for further information



Service and Support

Users of our systems can benefit from our comprehensive, fully inclusive service and support.

We provide complete service and support for all of our customers to give reassurance that if things go wrong or you need expert advice, help is only an e-mail or phone call away.



Validation Services

Our Validation Service enables you to implement and get maximum value from your investments as soon as possible.

We work as a partner with your Quality Manager, System Manager and users to provide a tailored Validation Plan suited to your needs. Our Validation Specialists incorporate years of experience in GLP system validation, detailed knowledge of our systems, together with other industry standard systems to help you meet your company's requirements.



Training

LabLogic can provide a variety of training courses and workshops to help you get the most out of your instrument and software.

All training is performed by our expert Product and Support Specialists who have many years experience in the development and use of the instruments and software.

Certificates can be provided to complement your internal GLP training records.

Related Products

Hidex Triathler

Liquid Scintillation Counter and Gamma Counter



Hidex 300 SL

Liquid Scintillation Counter



Posi-RAM™

PET Metabolite radio-HPLC Detector



Scan-RAM™

PET/SPECT radio-TLC Scanner



Flow-RAM™

PET/SPECT radio-HPLC Detector



Consumables



USA & Canada

LabLogic Systems, Inc.

East Pointe Park, 1040 East Brandon Blvd.
Brandon, FL 33511-5509, USA

E-mail: solutions@lablogic.com

Tel: +1-813-626-6848

Fax: +1-813-620-3708

Web: www.lablogic.com

Europe & Worldwide

LabLogic Systems Limited

Paradigm House, 3 Melbourne Avenue
Broomhill, Sheffield, S10 2QJ, UK

E-mail: solutions@lablogic.com

Tel: +44 (0)114 266 7267

Fax: +44 (0)114 266 3944

Web: www.lablogic.com



Certificate No: 10926
ISO 9001



Certificate No: 1535
ISO 9001



EXPERIENCE & EXPERTISE